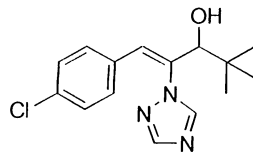
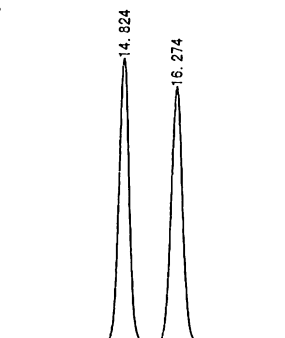
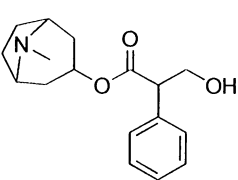
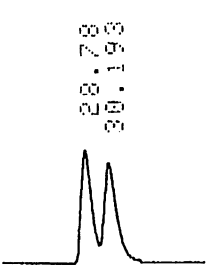
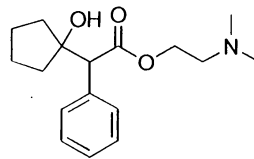
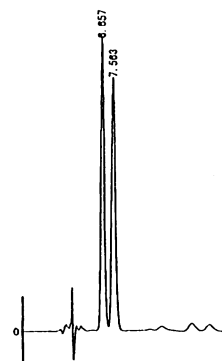
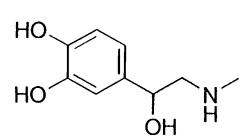
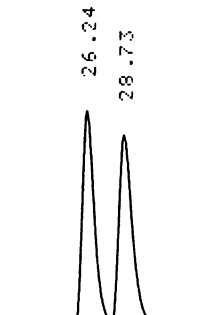
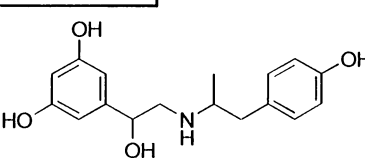
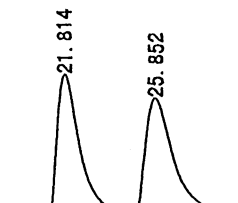
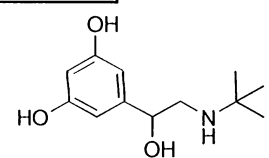
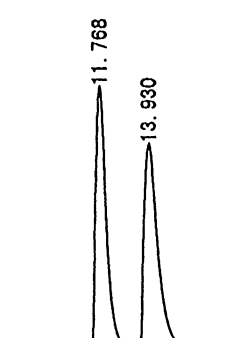
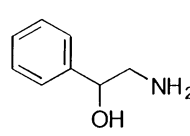
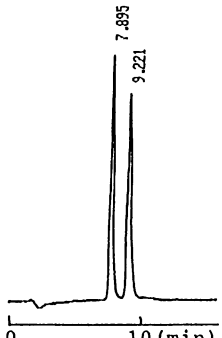


B-09	Uniconazole	B-10 Atropine
 OA-4700 Column:4.6mmI.D.×25cm Mobile phase:hexane/ethanol (98:2) Flow rate:1.0mL/min Detector:UV254nm Injection:1μL (1mg/mL)		 OA-4500 Column:4.6mmI.D.×25cm×2 Mobile phase:hexane/1,2-dichloroethane/ methanol/trifluoroacetic acid (250:140:10:1) Flow rate:1.0mL/min Detector:UV254nm
OA-4600 1.14 OA-4700 1.13 OA-4000 1.10 { R, S }	OA-4100 1.10 { R, S }	 OA-4500 1.06 AGP 1.29
B-11	Cyclopentolate	
 OA-7000 Column:4.6mmI.D.×25cm Mobile phase:40% acetonitrile in 20mmol/L pot.ph.b. pH3.0 Flow rate:0.7mL/min Detector:UV254nm Injection:3μL (3.5mg/mL)		
OA-7000 1.47 OA-7100 1.18 OA-4100 1.09	OA-4700 1.09 OA-4400 1.07 OA-4500 1.07	OA-4800 1.06 AGP 1.35
B-12	Epinephrine	B-13 Fenoterol
 OA-4900 Column:4.6mmI.D.×25cm Mobile phase:hexane/tetrahydrofuran/ trifluoroacetic acid (270:100:30:1) Flow rate:1.0mL/min Detector:UV254nm Injection:1μL (5mg/mL)		 OA-4900 Column:4.6mmI.D.×25cm Mobile phase:hexane/tetrahydrofuran/ methanol/trifluoroacetic acid (60:35:5:0.2) Flow rate:1.0mL/min Detector:UV280nm Injection:5μL (1mg/mL)
OA-4900 1.11 CBH 1.74		 OA-4900 1.22 OA-4500 1.13 OA-4700 1.10
B-14	Terbutaline	B-15 2-Amino-1-phenylethanol
 OA-4900 Column:4.6mmI.D.×25cm Mobile phase:hexane/tetrahydrofuran/ methanol/trifluoroacetic acid (60:35:5:0.2) Flow rate:1.0mL/min Detector:UV280nm Injection:5μL (1mg/mL)		 OA-5500 Column:4.6mmI.D.×15cm Mobile phase:1mmol/L copper (II) sulfate in water Flow rate:1.0mL/min Detector:UV254nm
OA-4900 1.25 OA-4700 1.13 AGP 1.31		 OA-5500 1.19 CBH 1.10