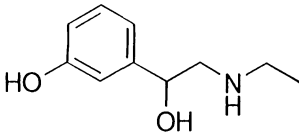
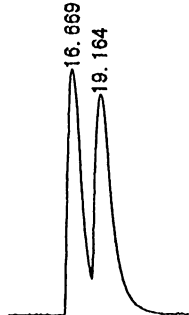
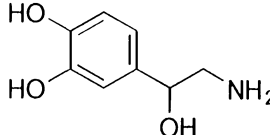
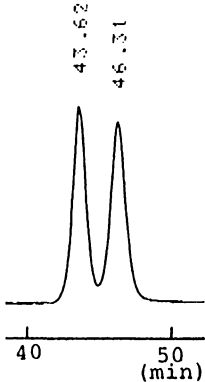
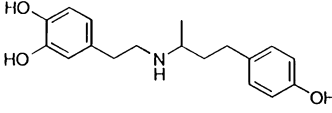
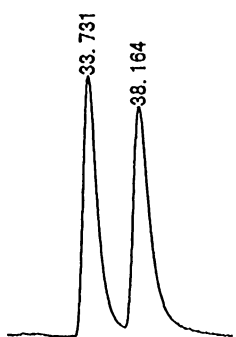
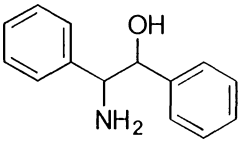
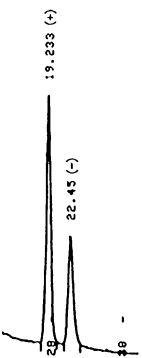
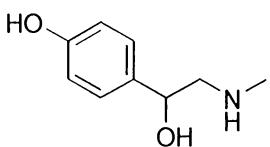
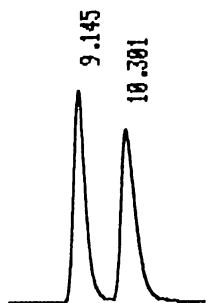
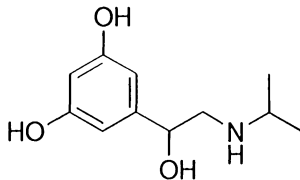
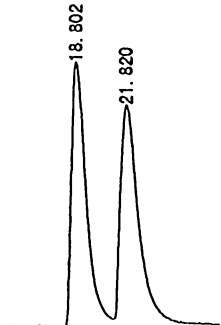
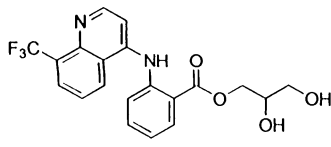
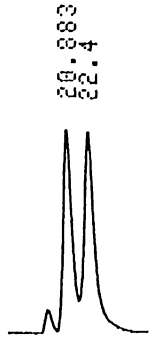
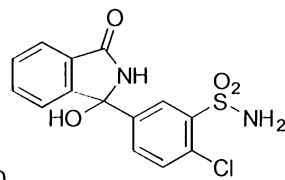
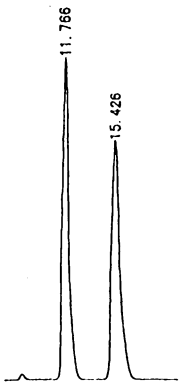


B-24 Etilefrine  OA-4900 Column:4.6mmI.D.×25cm Mobile phase:hexane/tetrahydrofuran/ methanol/trifluoroacetic acid (60:35:5:0.2) Flow rate:1.0mL/min Detector:UV280nm Injection:5μL (1mg/mL)  OA-4900 1.19	B-25 Norepinephrine  OA-4400 Column:4.6mmI.D.×25cm Mobile phase:hexane/tetrahydrofuran/ methanol/trifluoroacetic acid (270:100:30:1) Flow rate:1.0mL/min Detector:UV280nm Injection:1μL (5mg/mL)  OA-4400 1.07 CBH 2.13
B-26 Dobutamine  OA-4800 Column:4.6mmI.D.×25cm Mobile phase:hexane/tetrahydrofuran/ methanol/trifluoroacetic acid (60:35:5:0.2) Flow rate:1.0mL/min Detector:UV280nm Injection:2μL (1mg/mL)  OA-4800 1.15 CBH 1.81	B-27 erythro-2-Amino-1,2-diphenylethanol  OA-4900 Column:4.6mmI.D.×25cm Mobile phase:hexane/2-propanol/methanol/ trifluoroacetic acid (100:10:5:0.1) Flow rate:1.0mL/min Detector:UV254nm  OA-4900 1.20 { +, - } OA-4500 1.14 { +, - } OA-5500 1.13 OA-7000 1.10
B-28 Synephrine  OA-5000 Column:4.6mmI.D.×15cm Mobile phase:1mmol/L copper (II) acetate and 20mmol/L ammonium acetate in water (pH6.5) Flow rate:0.8mL/min Detector:UV254nm Injection:3μL (4mg/mL)  OA-5000 1.19 OA-4900 1.10	B-29 Metaproterenol  OA-4900 Column:4.0mmI.D.×25cm Mobile phase:hexane/tetrahydrofuran/ methanol/trifluoroacetic acid (60:35:5:0.2) Flow rate:1.0mL/min Detector:UV280nm Injection:5μL (1mg/mL)  OA-4900 1.19 OA-7000 1.12
B-31 Floctafenine  OA-4100 Column:4.6mmI.D.×25cm Mobile phase: hexane/tetrahydrofuran/ethanol (80:18:2) Flow rate:1.0mL/min Detector:UV254nm  OA-4100 1.09	B-32 Chlorthalidone  OA-7000 Column:4.6mmI.D.×25cm Mobile phase:20% acetonitrile in 20mmol/L pot.ph.b. pH3.0 Flow rate:0.7mL/min Detector:UV254nm  OA-7000 1.52 OA-7100 1.23 OA-4600 1.10 OA-4000 1.08 OA-4400 1.07 AGP 2.07