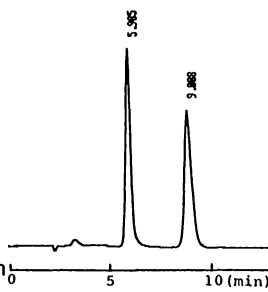
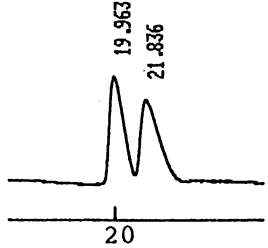
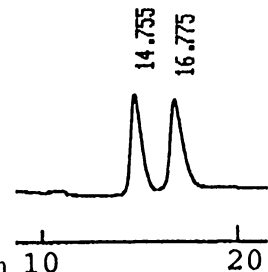
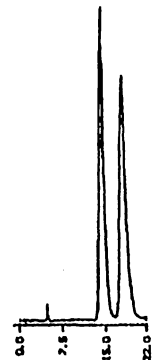
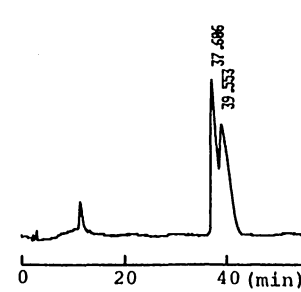
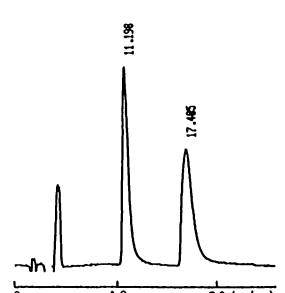
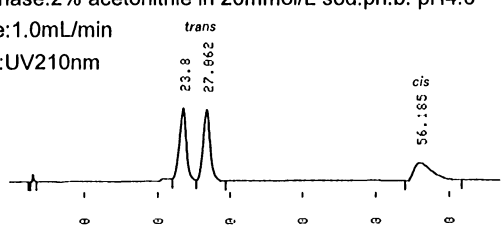
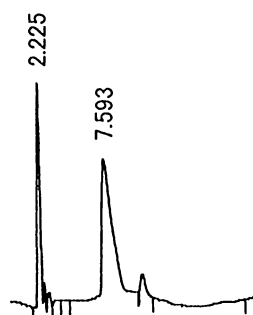


H-11 <chem>OC(=O)C1CCOC1</chem> OA-6000 Column:4.6mmI.D.×15cm Mobile phase:2mmol/L copper(II) sulfate in water/acetonitrile (90:10) Flow rate:1.0mL/min Detector:UV254nm Injection:3μL (6mg/mL)  OA-60001.70 OA-61001.28 OA-50001.06	H-12 <chem>OC(=O)C1CCOC1</chem> OA-6000 Column:4.6mmI.D.×15cm Mobile phase: 1mmol/L copper(II) sulfate in water Flow rate:1.0mL/min Detector:UV254nm Injection:3μL (8mg/mL)  OA-60001.10 OA-50001.05
H-13 <chem>OC(=O)C1C=CCOCC1</chem> OA-6000 Column:4.6mmI.D.×15cm Mobile phase:2mmol/L copper(II) sulfate in water/acetonitrile (90:10) Flow rate:1.0mL/min Detector:UV254nm Injection:5μL (4mg/mL)  OA-61001.18 OA-60001.15 OA-50001.09 OA-32001.07	H-14 <chem>CC(C)(C)C(=O)C=C(C)C(O)C1=CC(=O)C=C(C)C1</chem> AGP Column:4.0mmI.D.×10cm Mobile phase: 75mmol/L sod.ph.b. pH5.0 Flow rate:0.9mL/min Detector:UV225nm  OA-47001.06 AGP1.26
H-15 <chem>OC(=O)C(Br)C</chem> OA-5000 Column:4.6mmI.D.×15cm Mobile phase:2mmol/L copper(II) sulfate in water/acetonitrile (85:15) Flow rate:1.0mL/min Detector:UV254nm  OA-50001.05	H-16 <chem>OC(=O)C(O)C(C)(C)CO</chem> OA-5000 Column:4.6mmI.D.×15cm Mobile phase:2mmol/L copper(II) sulfate in water/2-propanol (90:10) Flow rate:1.0mL/min Detector:UV254nm Injection:5μL (8mg/mL)  OA-50001.63
H-17 <chem>OC(=O)C1CCCCC1C(=O)O</chem> OA-7500 Column:4.6mmI.D.×25cm Mobile phase:2% acetonitrile in 20mmol/L sod.ph.b. pH4.0 Flow rate:1.0mL/min Detector:UV210nm  OA-75001.16	H-18 <chem>OC(=O)C(Cl)C</chem> OA-6100 Column:4.6mmI.D.×15cm Mobile phase: 2mmol/L copper(II) sulfate in water Flow rate:1.0mL/min Detector:UV254nm  OA-610018.60