


# 食品関連化合物

## アクリルアミド (Acrylamide)



アクリルアミド  
(Acrylamide)

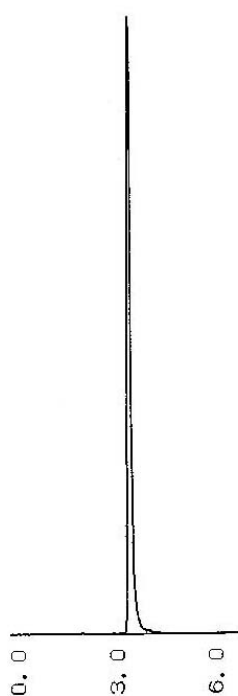

Column : SUMIPAX ODS D-210SLP ( 3  $\mu$  m, 2mm  $\times$  150mm )  
Mobile phase : water / methanol ( 99/1 )  
Flow rate : 0.2 mL/min  
Pressure : 11.2 MPa  
Temperature : 25  
Detector : UV 254nm  
Injection : 0.5  $\mu$  L ( 2.5mg/mL, methanol )