



その他

ヒドロキシ安息香酸 (Hydroxybenzoic acids)

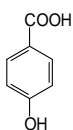

p-ヒドロキシ安息香酸
(*p*-Hydroxybenzoic acid)

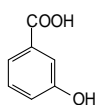


m-ヒドロキシ安息香酸
(*m*-Hydroxybenzoic acid)




安息香酸
(Benzoic acid)




サリチル酸
(Salicylic acid)

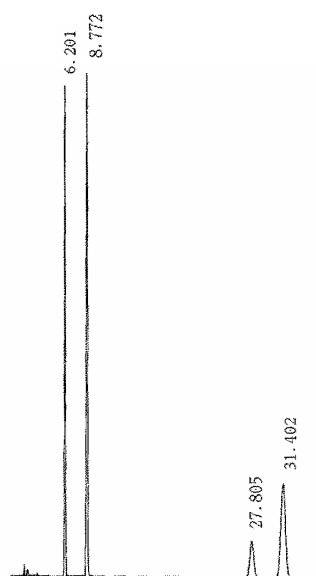

SUMIPAX ODS A-217(5 μm)

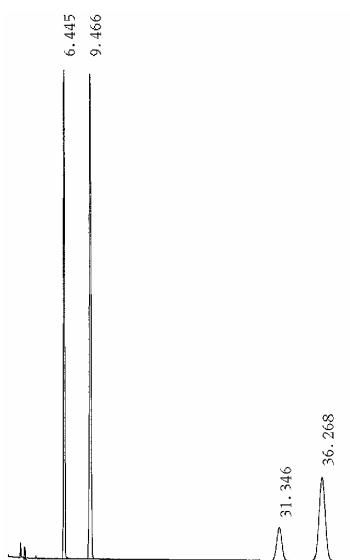

SUMIPAX ODS C-217(5 μm)

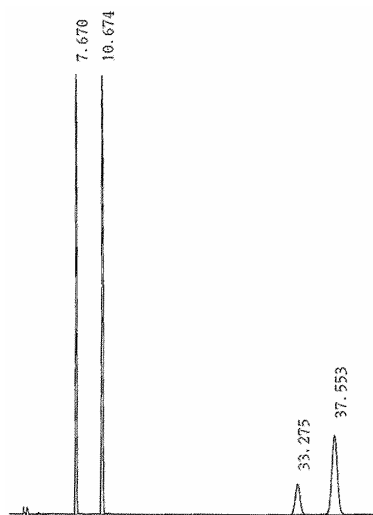

SUMIPAX ODS D-210SLP(3 μm)

Column : 4.6mm × 150mm
 Mobile phase : acetonitrile / water / acetic acid (5/90/5)
 Flow rate : 1.0 mL/min
 Temperature : 37
 Detector : UV 285nm
 Injection : 10 μL (0.1mg/mL respectively, mobile phase)